



LABCOLD™

medical & scientific refrigeration

Portable Vaccine Carrier from LABCOLD



The Labcold portable vaccine and sample carrier is designed to solve the problem of cold chain compliance on the move as it can be powered by both mains and vehicular power sources.

There is an easy to read digital temperature display and because it is accurate to one decimal place, the Labcold vaccine carrier can be calibrated so you can be confident that the temperature displayed is as experienced by the contents of the carrier.

Using a thermoelectric, solid state cooling system there are no dangerous liquids that could escape or moving parts that could get damaged.

- ❄ Supplied with two power leads. 230 volt for the mains and the other designed to fit 12v power sockets in vehicles
- ❄ Precise digital controller with real time and set point temperature display means the carrier can be validated in accordance with DH Green Book chapter 3 requirements
- ❄ Smooth, easy clean plastic exterior, aluminium interior
- ❄ Compact and lightweight

The carrier comes complete with a convenient shoulder strap and the interior features a divider and cushioned lids to protect the contents in transit.

Purpose designed for any professional who needs to give vaccines in the field, this carrier can be used in conjunction with Labcold data loggers for extra temperature security.


Please be aware this product is not designed to be run for extended periods of time nor to hold temperature when not connected to power source.

Not designed for the transit of blood products

Part Number	RPDF0012E
Temperature Range	+2°C to +8°C
Exterior Dimensions (HxWxD)	290 x 420 x 290mm
Interior Dimensions (HxWxD)	160 x 330 x 160mm
Weight	6.2kg
Energy Consumption	1.3 kWh/24hr
Carry Strap	✓
Digital temperature control	✓
Temperature display	✓
Moulded feet for stability in transit	✓
12 volt DC and 240 AC volt operation	✓



Handy compartment for vehicle cable

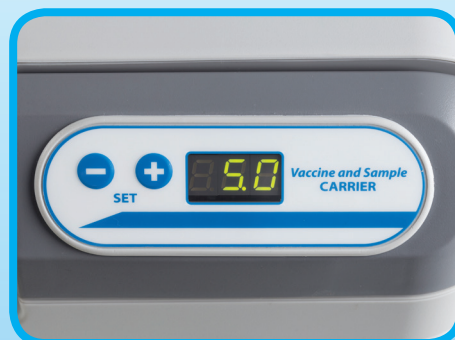

Clear digital temperature display



Hygienic easy clean exterior





Wolf Laboratories Limited

www.wolflabs.co.uk

Wolflabs

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

