



RLFF13246 Sparkfree Fridge Freezer



The **Labcold** Sparkfree range is engineered to store flammable materials safely with no sources of ignition inside the chambers. Anything that could cause a spark is located outside of the chamber.

Warning: Do not store open containers of volatile substances in this fridge freezer.

Did you know that 25% of chemical accidents are caused by improper storage?

Source: US Chemical Safety and Laboratory Investigation Board



Designed specifically for laboratory safety with a sparkfree interior and an optional factory fitted door lock for optimum security (ask for lock model RLFF13246-LK).

This fridge freezer provides a large storage capacity on a small floor area with two refrigeration systems so you can control the temperature of the fridge and freezer separately.

Part Number	RLFF13246
Capacity (Litres)	Fridge - 228
Capacity (Litres)	Freezer - 134
Style	Freestanding
Exterior Dimensions (HxWxD)	1965 x 595 x 630mm
Weight	80.4kg
Energy Consumption	1.36kWh/24hr
Operating Temperature	Fridge +2°C to +10°C Freezer -18°C to -23°C
Factory Set Point	Fridge +4°C Freezer -20°C
Interior Construction	Moulded Plastic
Exterior Construction	White Enamel
Door Style	Solid
Door Hinging	Right
Door Lock	Lock Model RLFF13246-LK
Shelves	Fridge - 4
Drawers	Freezer - 4
Feet	4
Temperature Display	Optional
High Temperature Alarm	Optional
Low Temperature Alarm	Optional
Defrost	Fridge - Auto Freezer - Manual
Eco Friendly Refrigerant	✓





Wolflabs


Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk

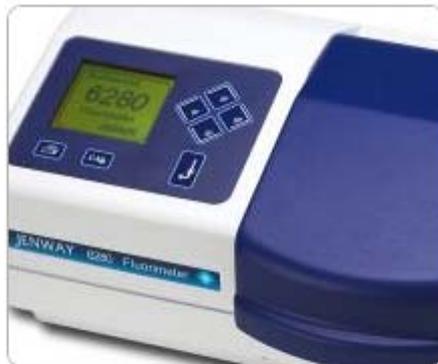

Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

