



# LABCOLD™

medical & scientific refrigeration



## Labcold -40 to -86°C Ultra Low Temperature 500 litre Freezer

Part Number	LULT0500U
Capacity (Litres)	500
Style	Freestanding
Exterior Dimensions (HxWxD)	1970 x 840 x 880mm
Interior Dimensions (HxWxD)	1310 x 600 x 630mm
Weight	320kg
Operating Temperature	-40°C to -86°C
Interior Construction	Stainless Steel
Exterior Construction	Epoxy Painted Steel
Door Style	Solid
Door Hinging	Right
Remote Alarm contacts	✓
High Temperature Alarm	✓
Low Temperature Alarm	✓
Mains Failure Alarm	✓
Low Battery Alarm	✓
Condenser High Temperature alarm	✓
Digital Temperature Display	✓
Key Operated Door Lock	✓
Password Protected Temperature Settings	✓
Battery Back up (Display and Alarms)	✓
Shelves	3
Inner Doors	4
Castors	4
Access Port	✓
Eco Friendly Refrigerant	✓
5 year warranty 3 x years parts and labour 2 x years parts only	✓

Labcold Ultra Low Temperature freezers use the tried and trusted cascade refrigeration and have been designed to occupy less than 1m<sup>2</sup> of floor space.



Complete with an easy to use digital controller conveniently situated on the front of the freezer, the temperature settings can be password protected and it always displays the set point for extra peace of mind. It also features alarm battery

back -up so should there be a mains power failure the freezer will alarm immediately to give you as much time as possible to deal with the emergency.

Backed by a 5 year parts and 3 year labour warranty, this Ultra Low Temperature freezer from Labcold is designed to be a robust and reliable place to store your samples.

We can also supply racking, why not ask for details at time of order?





**Wolflabs**


# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)

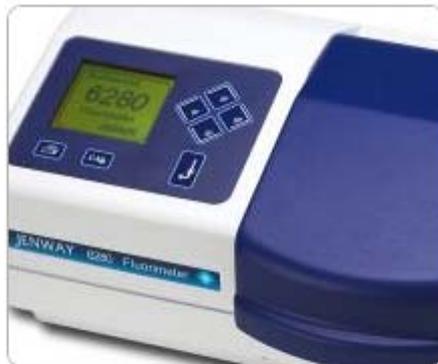

**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

