

The New Entris<sup>®</sup>.  
Just Weighing.




# The New Entris® . Just Weighing.

Receive the reliability and quality of a Sartorius balance, but just pay for the features you really need.

The new Entris® has been specially designed to achieve effective and reliable weighing results in your daily work. The easy and clearly structured user interface, with its logical key assignment and excellent readability, are ideal features for safe operation.


## Easy Cleaning

The top and side sliding doors and easy-to-remove stainless steel parts of the draft shield chamber, make Entris® convenient for thorough cleaning.

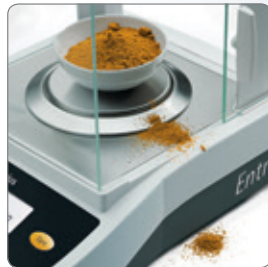

## Internal Adjustment

Each model is also available with an internal, motorized adjustment function for best weighing accuracy and operating convenience.

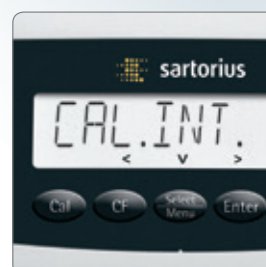



## High Accuracy

The monolithic weighing system ensures long-term high measurement accuracy and represents the proverbial reliability and quality of a Sartorius balance.



## Basic Applications

Entris® supports you in your daily work with integrated basic applications: animal weighing, density determination, counting, conversion or percentage.



## User-friendly Operation

Excellent readability due to the backlit and high-contrast display. The easy navigation with function keys and the simple-to-read level indicator on the front enable an user-friendly and fast operation.



## Technical Data Entris®

Model	Weighing capacity (g)	Readability (mg)	Repeatability (mg)	Linearity (mg)	Typical stabilization time (s)	Weighing pan size (mm)	Dimensions D×W×H (mm)	Calibration
<b>Analytical balances</b>								
ENTRIS224i-1S	220	0.1	0.1	0.2	2.5	∅ 90	230×303×330	Internal
ENTRIS224-1S	220	0.1	0.1	0.2	2.5	∅ 90	230×303×330	External
ENTRIS124i-1S	120	0.1	0.1	0.2	2.5	∅ 90	230×303×330	Internal
ENTRIS124-1S	120	0.1	0.1	0.2	2.5	∅ 90	230×303×330	External
ENTRIS64i-1S	60	0.1	0.1	0.2	2.5	∅ 90	230×303×330	Internal
ENTRIS64-1S	60	0.1	0.1	0.2	2.5	∅ 90	230×303×330	External
<b>Precision balances</b>								
ENTRIS623i-1S	620	1	1	2	1.0	∅ 115	230×303×136	Internal
ENTRIS623-1S	620	1	1	2	1.0	∅ 115	230×303×136	External
ENTRIS423i-1S	420	1	1	2	1.0	∅ 115	230×303×136	Internal
ENTRIS423-1S	420	1	1	2	1.0	∅ 115	230×303×136	External
ENTRIS323i-1S	320	1	1	2	1.1	∅ 115	230×303×136	Internal
ENTRIS323-1S	320	1	1	2	1.1	∅ 115	230×303×136	External
ENTRIS153i-1S	150	1	1	2	1.3	∅ 115	230×303×136	Internal
ENTRIS153-1S	150	1	1	2	1.3	∅ 115	230×303×136	External
ENTRIS6202i-1S	6,200	10	10	30	1.5	180×180	230×303×91	Internal
ENTRIS6202-1S	6,200	10	10	30	1.5	180×180	230×303×91	External
ENTRIS4202i-1S	4,200	10	10	30	1.5	180×180	230×303×91	Internal
ENTRIS4202-1S	4,200	10	10	30	1.5	180×180	230×303×91	External
ENTRIS3202i-1S	3,200	10	10	30	1.5	180×180	230×303×91	Internal
ENTRIS3202-1S	3,200	10	10	30	1.5	180×180	230×303×91	External
ENTRIS2202i-1S	2,200	10	10	30	1.5	180×180	230×303×91	Internal
ENTRIS2202-1S	2,200	10	10	30	1.5	180×180	230×303×91	External
ENTRIS822i-1S	820	10	10	30	1.5	∅ 150	230×303×87	Internal
ENTRIS822-1S	820	10	10	30	1.5	∅ 150	230×303×87	External
ENTRIS8201i-1S	8,200	100	100	300	1.5	180×180	230×303×91	Internal
ENTRIS8201-1S	8,200	100	100	300	1.5	180×180	230×303×91	External
ENTRIS5201i-1S	5,200	100	100	300	1.5	180×180	230×303×91	Internal
ENTRIS5201-1S	5,200	100	100	300	1.5	180×180	230×303×91	External
ENTRIS2201i-1S	2,200	100	100	300	1.5	180×180	230×303×91	Internal
ENTRIS2201-1S	2,200	100	100	300	1.5	180×180	230×303×91	External

For accessories see Entris® list of accessories.





**Wolflabs**

# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax:01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

